

Introducing The New Revolutionary Wireless STATPATCH Holter System

MOBILE HOLTER MONITORING CAN ADD UP TO **\$10,000**/YEAR PER UNIT IN REVENUE BASED ON 2 STUDIES PER WEEK

The StatPatch is a full 12-Lead ECG that is 100% bluetooth wireless and FDA - approved. **No Equipment Purchase Required, no out of pocket expense.** Quick and easy setup.

- ♦ Minimal preparation required
- ♦ System provides optimal patient comfort and acceptance
- ♦ Patient return not required

BILLING OPTIONS:


- Option 1.** No Charge to physician for the StatPatch Wireless System.
- Physician Bills CPT 93225 for hookup and recording and 93227 for review of report.
reimbursement totals \$79 (based on Florida 2007 Medicare rate)
- Option 2.** Physician purchases the StatPatch Wireless System.
- Physician Bills CPT 93230 which is the global 24 hour procedure.
Reimbursement totals \$165 (based on Florida 2007 Medicare rate)



Each patch is a single, self-contained, self-adhesive unit containing electrodes with pre-filled electrolyte gel, a processing board, and lithium battery. Minimal knowledge of anatomy is required for placement. Attachment is as easy as removing the backing from the self-adhesive StatPatch and aligning the patch with the top, touching the inter-clavicular space and the bottom resting between the fourth and fifth intercostal space. **No shaving of chest hair is necessary** (only in extreme cases). The ergonomic design of the STAT-PATCH™ provides optimal patient comfort without the worry of electrodes falling off.

Real-Time Data

Utilizing Bluetooth® technology, the STATPATCH™ transmits ECG data to a smartphone that patients either carry with them or place within 15 feet of themselves. Data is transmitted by the smartphone to an independent diagnostic testing facility that stores, and processes this data. Healthcare providers can rapidly access reports via a dedicated HIPAA (Health Insurance Portability and Accountability Act)-compliant Web site.

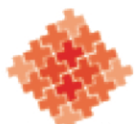

Telzuit

Medical Technologies, Inc.
www.Telzuit.com



For more Information or for a demonstration, call New Line Medical at **1-800-452-8909**
www.NewLineMedical.com