



The Effects of Irrigants on the Survival of Human Stem Cells of the Apical Papilla, Including Endocyn

Mark B. Scott, II, DMD,* Gregory S. Zilinski, DDS,*[†] Timothy C. Kirkpatrick, DDS,[†] Van T. Himel, DDS,* Kent A. Sabey, DDS,* and Thomas E. Lallier, PhD[‡]

Abstract

Introduction: Endocyn, a pH-neutral solution of hypochlorous acid and hypochlorite has been developed for use as an endodontic irrigant. The purpose of this study was to evaluate the effect of Endocyn on human periodontal ligament (PDL) fibroblasts, rat osteosarcoma cells (UMR-106), and stem cells of the apical papilla (SCAP) compared with other commonly used endodontic irrigants. **Methods:** To determine cytotoxicity, cells were exposed to various concentrations of Endocyn, 6% sodium hypochlorite (NaOCl), 17% EDTA, and 2% chlorhexidine for 10 minutes, 1 hour, or 24 hours. Cell survival was measured fluorescently using calcein AM. Endocyn also was tested for its ability to inhibit SCAP proliferation and alkaline phosphatase activity. Finally, SCAP transcript expression was examined via reverse-transcriptase polymerase chain reaction. **Results:** Endocyn was no more toxic to PDL and UMR cells than water for up to 24 hours. Endocyn concentrations of 50% were toxic to SCAP after 1 hour of exposure. Endocyn concentrations of >20% inhibited SCAP proliferation, whereas concentrations of ≥10% inhibited alkaline phosphatase activity. Exposure of SCAP to 10% Endocyn for 3 days did not alter most transcript expression, but did significantly reduce the expression of alkaline phosphatase, fibromodulin, and osteomodulin. **Conclusion:** Endocyn was significantly less cytotoxic to PDL, UMR-106, and SCAP cells compared with other commonly used endodontic irrigants. High concentrations of Endocyn did inhibit some transcript expression and alkaline phosphatase activity, indicating a potential reduction in the osteogenic potential of stems cells exposed to Endocyn. (*J Endod* 2018;44:263–268)

Key Words

Cell survival, chlorhexidine, hypochlorous acid, stem cells, sodium hypochlorite

It is well established that microorganisms are the principal cause of endodontic disease (1–3), thus the paramount goal of endodontic treatment is their elimination.

Regenerative endodontics

is indicated for management of immature permanent teeth with pulpal necrosis. This procedure is performed with the goal of allowing and promoting further root development. This is accomplished by first disinfecting the canal space, and then allowing the ingress of stem cells into that space. Those stems cells of the apical papilla (SCAP) typically reside in close proximity to the root end and have exhibited faster proliferation, and greater differentiation into odontoblast-like cells than other dental stem cells (4, 5). SCAP survival during the disinfection process is essential to the regenerative process (6). If these stem cells are exposed to standard irrigants or medicaments, it may render them nonviable, less functional, or unable to migrate appropriately to the needed site (7, 8). This creates a conundrum regarding one's ability to adequately disinfect the canal space while promoting the survival of native cells, allowing regeneration to proceed. Sodium hypochlorite (NaOCl) has proven to be an effective irrigant because of its ability to effectively disrupt microbial biofilms (9) and dissolve necrotic tissue (10). However, a significant disadvantage of using NaOCl is its toxicity to periodontal ligament (PDL) cells (11) and SCAPs, (7) and the significant pain and morbidity involved when it is extruded beyond the confines of the tooth (12). Chlorhexidine (CHX) is another endodontic irrigant that efficiently eliminates microorganisms but also negatively impacts PDL (11) and gingival fibroblasts (13). An optimal irrigant for regenerative procedures would disinfect the canal space while not harming human cells in the vicinity (6). It has been shown that 17% EDTA promotes SCAP survival and attachment. In contrast, 2% CHX is detrimental to SCAP survival, and 6% NaOCl + EDTA decreased cell viability when compared with EDTA alone (7). Unfortunately, the use of EDTA alone is not sufficient to consistently produce a disinfected canal (14). Other intracanal medicaments, such as triple- and double-antibiotic pastes, were both detrimental to SCAP survival, whereas Ca(OH)₂ promoted SCAP survival and proliferation (8, 15). Thus, many root canal medicaments have been demonstrated to be detrimental to SCAP and are therefore not suitable for regenerative endodontics. The ideal endodontic irrigant would eliminate microorganisms but have insignificant or no toxicity to healthy host cells (16).

Significance

Regenerative pulp therapy requires host stem cell survival in the root canal. Most endodontic irrigants are highly cytotoxic. Endocyn is less cytotoxic to stem cells of the apical papilla and is better suited for irrigation of teeth with open apices.

From the *Departments of Endodontics and [†]Cell Biology and Anatomy, Louisiana State University, School of Dentistry, New Orleans, Louisiana; and [‡]Department of Endodontics, Keesler Air Force Base, Biloxi, Mississippi.

Address requests for reprints to Thomas E. Lallier, Louisiana State University, Department of Cell Biology and Anatomy, 1100 Florida Avenue, New Orleans, LA 70119. E-mail address: tlalli@lsuhsc.edu
0099-2399/\$ - see front matter

Copyright © 2017 American Association of Endodontists.
<https://doi.org/10.1016/j.joen.2017.09.001>

Endocyn is a super-oxidized solution, manufactured from pure water and sodium chloride, shown to have antimicrobial properties (17, 18). According to the manufacturer, through a proprietary process, water and saltwater molecules are dissociated forming hypochlorous acid and hypochlorite ions (19). Microcyn, which is chemically similar to Endocyn, has been used as an antimicrobial wound treatment and surface disinfectant with a wide antimicrobial range in a hospital setting (18). Microcyn also has been shown to have antimicrobial properties, while also being nontoxic to human dermal fibroblast cultures when compared with hydrogen peroxide (20). It is also reported to be an effective treatment for mildly infected diabetic foot ulcers, with minimal adverse effects (21). Both Microcyn and Endocyn contain water, hypochlorous acid, and sodium chloride. However, Microcyn also contains NaOCl, whereas Endocyn also contains sodium sulfate and monobasic sodium phosphate. Thus, these solutions are similar, but not identical. Because Microcyn has shown to be a useful disinfectant while exhibiting minimal effects on human cells, Endocyn could potentially become a beneficial treatment in pulpal regeneration. The purpose of this study was to compare cell survival of PDL, osteosarcoma (UMR), and SCAP cells following exposure to Endocyn versus other traditional endodontic irrigants. This study also examined the effect of Endocyn on SCAP proliferation and gene expression.

Material and Methods

The SCAPs were obtained from Dr. Anibal Diogenes, Department of Endodontics, University of Texas Health Science Center at San Antonio. This SCAP cell line, RP-89, was derived from the apical papillae of a mandibular third molar from a single donor. RP-89 has shown to preferentially express and maintain mesenchymal stem cell markers in culture, coexpressing CD73, CD90, and CD105 in all passages evaluated (22). Expression of these molecular markers is considered the minimal criteria for the identification of true mesenchymal stem cells (23). Cells were maintained in minimum essential medium alpha (MEM α) containing 10% fetal bovine serum (FBS) at 37°C and 5% CO₂ (Invitrogen, GIBCO, Grand Island, NY).

Human PDL fibroblasts were obtained from patients with healthy gingiva and no periodontal disease who underwent extraction at the Department of Oral Surgery at the Louisiana State University School of Dentistry, New Orleans, LA (24). All tissues were obtained from subjects after informed written consent as prescribed in an approved institutional review board protocol (#004826). Cells were maintained in MEM α containing 10% FBS at 37°C and 5% CO₂. Cell lines of rat osteosarcoma cells (UMR-106) were obtained from American Type Culture Collection (ATCC-CRL-1661; Manassas, VA) and were maintained in Dulbecco's Modified Eagle's Medium containing 10% fetal calf serum and 200 U/mL penicillin and 200 mg/mL streptomycin (Invitrogen, GIBCO). In all cases, cells were removed from polystyrene substrates using 0.25% trypsin before use in subsequent experiments. Cells were centrifuged at 200g for 10 minutes, resuspended in media containing 10% serum, and cell numbers verified by hemocytometer.

Irrigants

Endocyn was obtained from Sonoma Pharmaceuticals, Inc (formerly Oculus Innovative Sciences), Petaluma, CA. Sodium hypochlorite (6% solution) and CHX (Vista-CHX Plus, 480) were purchased from Inter-Med, Inc/Vista Dental Products, Racine, WI. EDTA (17% solution) was purchased from Coltene/Whaledent, Inc, Cuyahoga Falls, OH.

Exposure to Irrigants

Cells were allowed to adhere for 24 hours to 48-well plates in MEM α containing 10% FBS. The cell media was then replaced with fresh MEM α containing various dilutions of endodontic irrigant obtained by diluting 6% NaOCl, distilled water, 2% CHX, 17% EDTA, and Endocyn. MEM α alone served as a negative control. Following exposure to the solutions for 10 minutes, 1 hour, or 24 hours, the cells were treated with calcein AM for 1 hour to render live cells fluorescent. Cells were rinsed twice with phosphate-buffered saline (PBS) and cell viability was determined using a Synergy 2 plate reader (Biotek Instruments, Inc, Winooski, VT), calibrated to measure fluorescence intensity with filters appropriate for approximately 480-nm excitation and approximately 520-nm emission. Values represent the means and standard deviations (SDs) of 8 samples (after subtraction of autofluorescence).

Cell Proliferation

SCAP cells were plated into 48-well tissue culture plate at 4000 cells in 200 μ L in MEM α , and allowed to adhere for 1 hour. Cells were incubated for up to 9 days at 37°C and adherent cells were rinsed with PBS and quantified fluorometrically using CyQuant fluorescent dye (Molecular Probes, Eugene, OR) (25). Fluorescence was measured using a Synergy 2 plate reader (Biotek Instruments, Inc) with filters appropriate for approximately 480-nm excitation and approximately 520-nm emission. Values represent the means and SDs of 8 samples (after subtraction of autofluorescence).

Alkaline Phosphatase Activity

SCAP cells exposed to Endocyn were assayed for alkaline phosphatase activity (25). Cells were plated into 48-well plates and exposed to various concentrations of Endocyn for up to 9 days. Cells were rinsed 3 times in PBS and frozen to lyse cells. P-nitrophenol phosphate 1 mg/mL in 0.1 M diethanolamine (pH 8.3) was added to each well (200 μ L per well) and incubated at 25°C for 30 minutes with gentle agitation. The enzymatic color reaction was stopped by the addition of 500 μ L 0.75 N NaOH, and assayed for 405-nm absorbance in a Biotek Synergy 2 plate reader (Biotek Instruments, Inc). Values represent the means and SDs of 8 samples.

Comparison of Transcript Expression

Differential transcript expression was examined by reverse-transcriptase polymerase chain reaction (PCR). Briefly, RNA was extracted from tissue and cell samples using Qiagen kits (Germantown, MD) and gene expression determined by reverse-transcriptase PCR. The number of cycles was chosen by creating a standard curve using serial dilutions of DNA template known to include the genes in question, and determining the PCR conditions that reveal twofold differences in gene expression over 3 orders of magnitude. Twenty-five to 35 cycles of PCR were chosen to be within the linear range of detection for all of the genes examined. The primers for these studies were derived from the published DNA sequences of human genes as described previously (24, 26).

Statistical Analysis

For quantitative analysis, 8 samples were averaged and the means and SDs were compared with the control values for untreated cells. Comparisons of experimental groups were performed via analysis of variance where normality of the data was confirmed using the Shapiro-Wilk normality test, with $P < .05$ considered to be significant.

The null hypothesis is that Endocyn is as cytotoxic as NaOCl or other endodontic irrigants.

Results

As anticipated, NaOCl and CHX were toxic to both PDL and UMR cells (Fig. 1), with a loss of 80% viability after exposure for only 10 minutes (Fig. 1A and B), and virtually no cell survival after 1 hour of exposure (Fig. 1C and D) at concentrations as low as 1%. In contrast, distilled water did not alter cell viability up to concentrations of 50%. Interestingly, Endocyn displayed no significant toxicity for these cells up to 50% ($P > .05$). In contrast, EDTA solutions of 1% and 5% displayed intermediate toxicity to these cells when exposed for up to 1 hour (Fig. 1C and D). At 24 hours, cell viability was again reduced

to nearly 0% (Fig. 1E and F). Thus, Endocyn was no more toxic to these cells than was water.

In addition, Endocyn promoted better SCAP cell survival than traditional irrigants (Fig. 2). Endocyn, at concentrations of $\leq 20\%$, was not significantly toxic to SCAP for all time periods ($P > .05$). Endocyn at 50% concentration showed slight but significant toxicity ($P < .05$) at 10 minutes with 80% SCAP survival (Fig. 2A), 20% SCAP survival at 1 hour (Fig. 2B), and 0% SCAP survival at 24 hours (Fig. 2C). Thus, Endocyn was significantly less cytotoxic to SCAP cells than all of the other endodontic irrigants tested, especially after a prolonged exposure of 24 hours.

Exposure of SCAP cells to Endocyn only significantly inhibited their proliferation when the concentration was 50% Endocyn (Fig. 3A). Similarly, high doses of Endocyn inhibited SCAP cell alkaline phosphatase activity (Fig. 3B). However, this became detectable only at

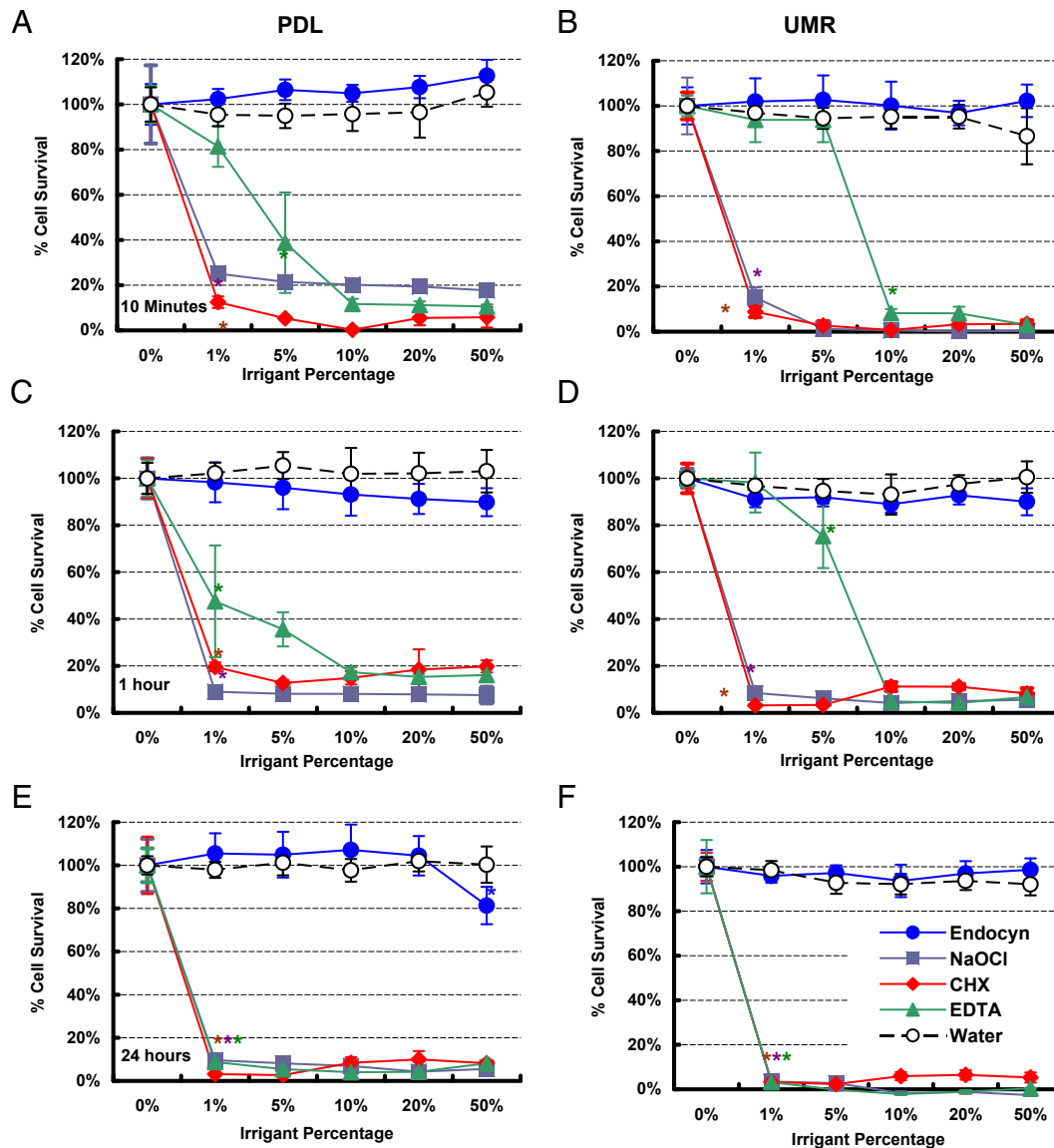


Figure 1. PDL and UMR cell survival is greater with Endocyn than other endodontic irrigants. PDL (A, C, and E) and UMR cells (B, D, and F) were plated for 24 hours before exposure to various dilutions of Endocyn, hypochlorite, EDTA, CHX, and distilled water. Live, calcein AM-labeled cells were examined after continuous treatment for 10 minutes (A and B), 1 hour (C and D), and 24 hours (E and F). Cell survival was measured fluorometrically, and the means and standard deviations of 8 samples were compared. Significant differences from untreated control cell survival are indicated by asterisks, $P < .05$. CHX, chlorhexidine; NaOCl, sodium hypochlorite; PDL, periodontal ligament; UMR, osteosarcoma cells.

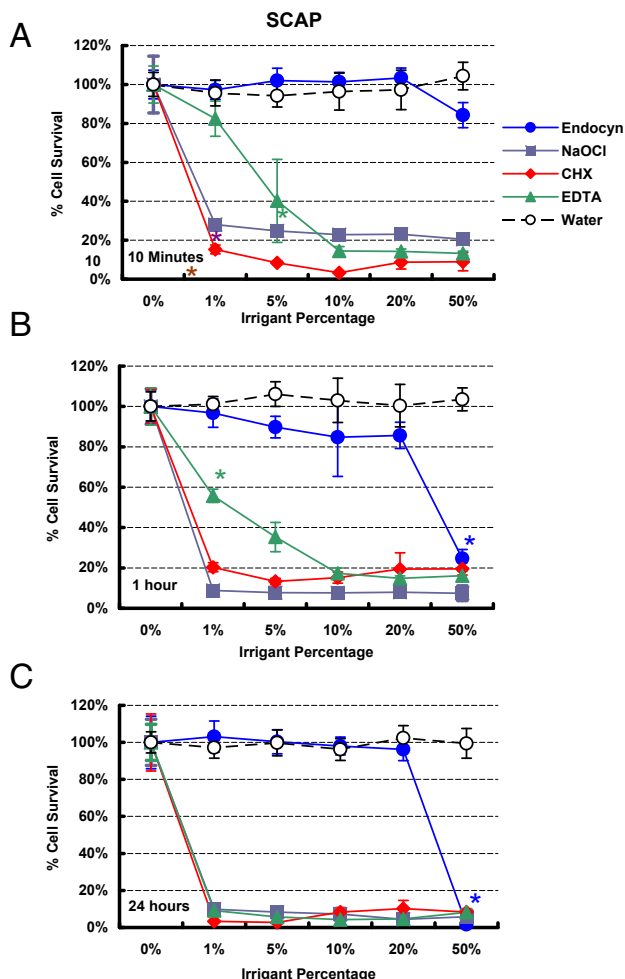


Figure 2. SCAP cell survival is greater with Endocyn than other endodontic irrigants. SCAP cells were plated for 24 hours before exposure to various dilutions of Endocyn, hypochlorite, EDTA, CHX, and distilled water. Live, calcein AM-labeled cells were examined after continuous treatment for 10 minutes (A), 1 hour (B), and 24 hours (C). Cell survival was measured fluorometrically, and the means and standard deviations of 8 samples were compared. Significant differences from untreated control cell survival are indicated by asterisks, $P < .05$. CHX, chlorhexidine; NaOCl, sodium hypochlorite; SCAP, stem cells of the apical papilla.

concentrations $\geq 10\%$ and after 7 days of exposure. Thus, only high doses of Endocyn appear to have detrimental effects on SCAP cells. Furthermore, an examination of SCAP transcript expression (Fig. 4) revealed that exposure with 10% Endocyn for 3 days significantly reduced alkaline phosphatase and 2 small leucine-rich proteoglycans (fibromodulin and osteomodulin) transcript expression. In contrast, Endocyn exposure has no appreciable effect on collagen Type I, integrin $\beta 1$, periostin, osteoprotegerin, osteonidogen, transforming growth factor $\beta 1$, growth differentiation factor 5, bone morphogenetic protein 2, bone morphogenetic protein 4, fibroblast growth factor 2, fibroblast growth factor 5, connective tissue growth factor, or matrix metalloproteinase 2. Thus, prolonged exposure to Endocyn selectively suppressed the expression of only a minority of transcripts examined.

Discussion

In this study, Endocyn demonstrated significantly less cytotoxicity to PDL, UMR-106, and SCAP cells compared with traditional endodontic

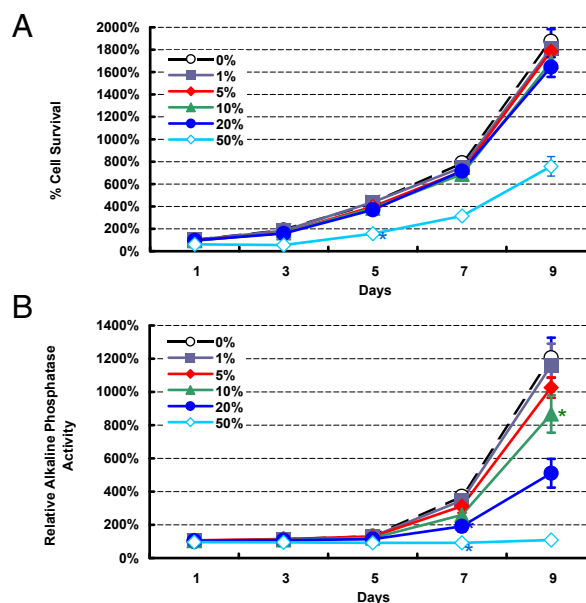


Figure 3. Endocyn inhibited SCAP cell proliferation and alkaline phosphatase activity at very high concentrations. SCAP cells were plated for 24 hours before exposure to various dilutions of Endocyn for up to 9 days. Cell proliferation (A) was measured fluorometrically using Cyquant and alkaline phosphatase activity was measured colorimetrically using P-nitrophenol phosphate (B). The means and standard deviations of 8 samples were compared. Significant differences from untreated control cell survival are indicated by asterisks, $P < .05$. SCAP, stem cells of the apical papilla.

irrigants. Endocyn reduced SCAP proliferation, alkaline phosphatase activity, and expression of fibromodulin, and osteomodulin transcripts. Fibromodulin and osteomodulin are small leucine-rich proteoglycans expressed by dental pulp cells (27) and PDL cells (24) and are involved in the fine architecture of the extracellular matrix (28, 29). However, Endocyn did not alter the expression of most transcripts examined, and therefore should be considered much less detrimental to these stem cells than traditional irrigants.

Ideal properties for endodontic irrigants include antimicrobial activity, the ability to dissolve necrotic tissue, to aid in the debridement of the canal system, and be nontoxic to the periradicular tissues (30). An ideal irrigant for regenerative procedures would disinfect the canal space while not harming human cells in the vicinity (16, 31). In this study, Endocyn proved to be less cytotoxic to SCAP than other tested irrigants and could prove useful not just for the disinfection of immature permanent teeth during regenerative procedures, but other endodontic procedures as well.

The active ingredient in Endocyn is hypochlorous acid, which has significant bactericidal activity due to its ability to penetrate bacterial cell membranes resulting in protein degradation (32). If Endocyn demonstrates the ability to eliminate microorganisms and disrupt microbial biofilms, it may become an additional option for irrigation during root canal therapy. Because of its relative nontoxicity, at low concentrations, Endocyn may become a safe irrigant option in cases whereby extrusion of irrigants is a concern, as in perforations or open apex cases. It could also prove useful in the dental education setting, as a safe alternative irrigant for clinicians with limited endodontic experience.

This study examined only the percent cell survival and proliferation *in vitro* for relatively short times (up to 9 days). Endocyn and its long-term effect on SCAP would also need to be studied. Cell proliferation, alkaline phosphatase activity, and expression of dentin-specific

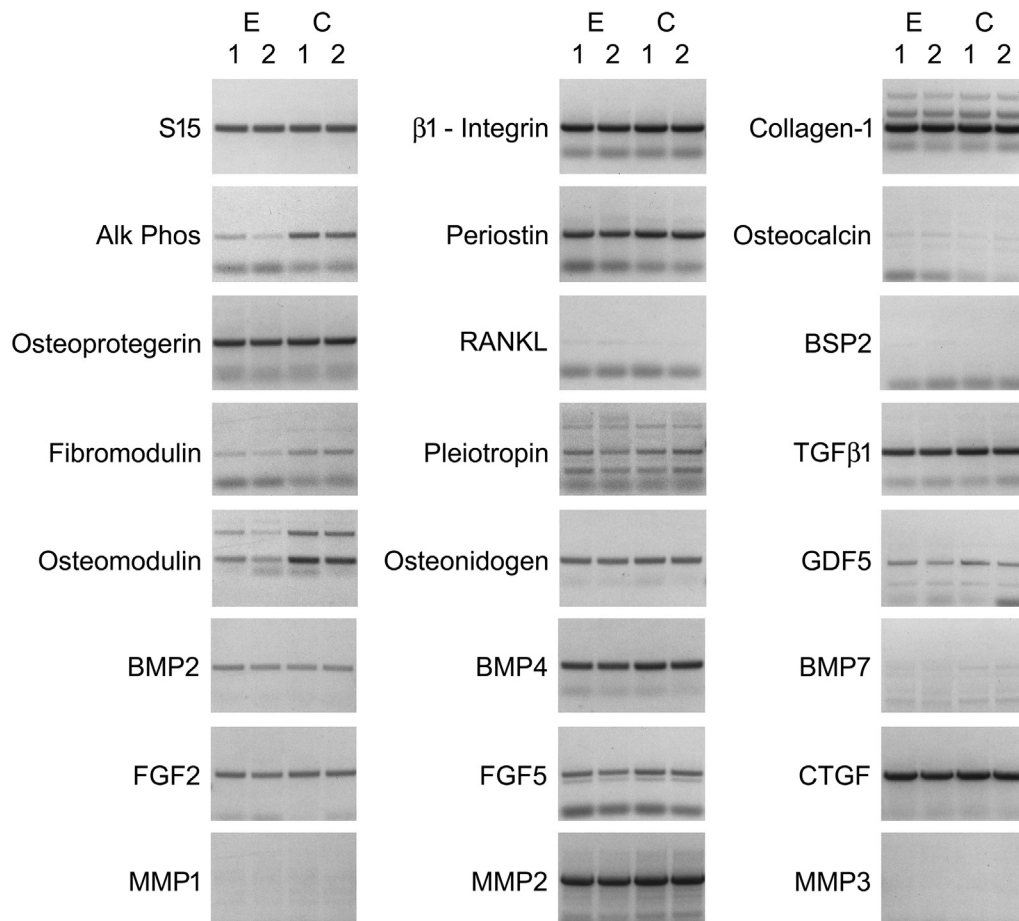


Figure 4. Endocyn selectively reduced SCAP expression of several osteogenic transcripts. SCAP cells were allowed to adhere for 24 hours before exposure to 10% Endocyn (E) or control growth media (C). After 3 days of exposure, transcript expression was assayed using gene-specific primers and RT-PCR. Only the expression of alkaline phosphatase, osteomodulin, and fibromodulin were inhibited by exposure to Endocyn. The ribosomal subunit (S15) was used as control for template integrity. Alk Phos, alkaline phosphatase; BMP, bone morphogenetic protein; BSP, bone sialoprotein; CTGF, connective tissue growth factor; FGF, fibroblast growth factor; GDF, growth differentiation factor; MMP, matrix metalloproteinase; RANKL, receptor activator of nuclear factor kappa-B ligand; RT-PCR, reverse-transcriptase polymer chain reaction; TGF, transforming growth factor; SCAP, stem cells of the apical papilla.

proteins would show SCAP survival and also whether they would continue to function properly.

Conclusion

Endocyn demonstrated less cytotoxicity to PDL, UMR-106, and SCAP cells than other traditional endodontic irrigants. Further testing regarding its ability to eliminate microorganisms, disrupt biofilms, and dissolve tissue, and its interactions with current endodontic irrigants is necessary. With more research and continued demonstration of ideal or beneficial irrigant properties, Endocyn, could become a safe and effective alternative to traditional endodontic irrigation, especially for endodontic regenerative therapies as well as during the training of undergraduate dental students.

References

1. Kakehashi S, Stanley HR, Fitzgerald RJ. The effects of surgical exposures of dental pulps in germ-free and conventional laboratory rats. *Oral Surg Oral Med Oral Pathol* 1965;20:340–9.
2. Moller AJ, Fabricius L, Dahlen G, et al. Influence on periapical tissues of indigenous oral bacteria and necrotic pulp tissue in monkeys. *Scand J Dent Res* 1981;89:475–84.
3. Fabricius L, Dahlen G, Sundqvist G, et al. Influence of residual bacteria on periapical tissue healing after chemomechanical treatment and root filling of experimentally infected monkey teeth. *Eur J Oral Sci* 2006;114:278–85.
4. Sonoyama W, Liu Y, Yamaza T, et al. Characterization of the apical papilla and its residing stem cells from human immature permanent teeth: a pilot study. *J Endod* 2008;34:166–71.
5. Lovelace TW, Henry MA, Hargreaves KM, et al. Evaluation of the delivery of mesenchymal stem cells into the root canal space of necrotic immature teeth after clinical regenerative endodontic procedure. *J Endod* 2011;37:133–8.
6. Huang GT, Sonoyama W, Liu Y, et al. The hidden treasure in apical papilla: the potential role in pulp/dentin regeneration and bioroot engineering. *J Endod* 2008;34:645–51.
7. Trevino EG, Patwardhan AN, Henry MA, et al. Effect of irrigants on the survival of human stem cells of the apical papilla in a platelet-rich plasma scaffold in human root tips. *J Endod* 2011;37:1109–15.
8. Ruparel NB, Teixeira FB, Ferraz CC, et al. Direct effect of intracanal medicaments on survival of stem cells of the apical papilla. *J Endod* 2012;38:1372–5.
9. Clegg MS, Vertucci FJ, Walker C, et al. The effect of exposure to irrigant solutions on apical dentin biofilms *in vitro*. *J Endod* 2006;32:434–7.
10. Hand RE, Smith ML, Harrison JW. Analysis of the effect of dilution on the necrotic tissue dissolution property of sodium hypochlorite. *J Endod* 1978;4:60–4.
11. Chang YC, Huang FM, Tai KW, et al. The effect of sodium hypochlorite and chlorhexidine on cultured human periodontal ligament cells. *Oral Surg Oral Med Oral Pathol Oral Radiol Endod* 2001;92:446–50.
12. Zhu WC, Gyamfi J, Niu LN, et al. Anatomy of sodium hypochlorite accidents involving facial ecchymosis—a review. *J Dent* 2013;41:935–48.

13. Tsourounakis I, Palaologou-Gallis AA, Stoute D, et al. Effect of essential oil and chlorhexidine mouthwashes on gingival fibroblast survival and migration. *J Periodontol* 2013;84:1211–20.
14. Wang Z, Shen Y, Haapasalo M. Effect of smear layer against disinfection protocols on *Enterococcus faecalis*-infected dentin. *J Endod* 2013;39:1395–400.
15. Althumairy RI, Teixeira FB, Diogenes A. Effect of dentin conditioning with intracanal medicaments on survival of stem cells of apical papilla. *J Endod* 2014;40:521–5.
16. Diogenes A, Henry MA, Teixeira FB, et al. An update on clinical regenerative endodontics. *Endod Topics* 2013;28:2–23.
17. Tanaka H, Hirakata Y, Kaku M, et al. Antimicrobial activity of superoxidized water. *J Hosp Infect* 1996;34:43–9.
18. Landa-Solis C, Gonzalez-Espinosa D, Guzman-Soriano B, et al. Microcyn: a novel super-oxidized water with neutral pH and disinfectant activity. *J Hosp Infect* 2005;61:291–9.
19. Oculus. *EndocynTM A New Product From the Family of MicrocynR*. California: Oculus Innovative Sciences; 2016.
20. Gonzalez-Espinosa D, Perez-Romano L, Guzman-Soriano B, et al. Effects of pH-neutral, super-oxidised solution on human dermal fibroblasts *in vitro*. *Int Wound J* 2007;4:241–50.
21. Landsman A, Blume PA, Jordan DA Jr, et al. An open-label, three-arm pilot study of the safety and efficacy of topical Microcyn Rx wound care versus oral levofloxacin versus combined therapy for mild diabetic foot infections. *J Am Podiatr Med Assoc* 2011;101:484–96.
22. Ruparel NB, de Almeida JF, Henry MA, et al. Characterization of a stem cell of apical papilla cell line: effect of passage on cellular phenotype. *J Endod* 2013;39:357–63.
23. Dominici M, Le Blanc K, Mueller I, et al. Minimal criteria for defining multipotent mesenchymal stromal cells. The International Society for Cellular Therapy position statement. *Cytotherapy* 2006;8:315–7.
24. Lallier TE, Spencer A, Fowler MM. Transcript profiling of periodontal fibroblasts and osteoblasts. *J Periodontol* 2005;76:1044–55.
25. Bonson S, Jeansonne BG, Lallier TE. Root-end filling materials alter fibroblast differentiation. *J Dent Res* 2004;83:408–13.
26. Bulmansi Z, Brady M, Stoute D, et al. Cigarette smoke extract induces select matrix metalloproteinases and integrin expression in periodontal ligament fibroblasts. *J Periodontol* 2012;83:787–96.
27. Buchaille R, Couble ML, Magloire H, et al. Expression of the small leucine-rich proteoglycan osteoadherin/osteomodulin in human dental pulp and developing rat teeth. *Bone* 2000;27:265–70.
28. Jan AT, Lee EJ, Choi I. Fibromodulin: a regulatory molecule maintaining cellular architecture for normal cellular function. *Int J Biochem Cell Biol* 2016;80:66–70.
29. Tashima T, Nagatoishi S, Sagara H, et al. Osteomodulin regulates diameter and alters shape of collagen fibrils. *Biochem Biophys Res Commun* 2015;463:292–6.
30. Harrison JW. Irrigation of the root canal system. *Dent Clin North Am* 1984;28:797–808.
31. Lin LM, Ricucci D, Huang GT. Regeneration of the dentine-pulp complex with revitalization/revascularization therapy: challenges and hopes. *Int Endod J* 2014;47:713–24.
32. Winter J, Ilbert M, Graf PC, et al. Bleach activates a redox-regulated chaperone by oxidative protein unfolding. *Cell* 2008;135:691–701.